

**T**he BagMixer® 400 line, the world reference in lab blenders, is now completed by 100 and 3500 ml models. The new BagMixer® range is the most complete and modern lab blender system on the market. To extract bacteria from a solid sample, just place the sample and liquid in a sterile BagFilter®, and then slide the bag into the BagMixer®. Blend for 30 to 60 seconds (for most samples). The juice containing the filtered extract may be pipeted and analysed without any risk of alteration or cross contamination. BagMixer® is ready for another sample in a new sterile bag. **This is the safest and most economical method for solid sample preparation.**

3 year warranty\*



## Main features

stainless steel body & blending chamber  
easy cleaning  
adjustable paddles  
silent

## Options

W (window door) or  
P (plain door)  
CoolMix®, WarmMix® :  
temp. regulated.

## BagMixer® 100 “MiniMix®”

Food, clinical/medical, chemical,  
agricultural, pharmacy, cosmetic...

Timer : 10 - 360 sec. or ∞, 6 - 9 strokes/sec.  
35 x 20 x 17.5 cm, 8 kg - 110 V / 220 V  
Capacity : 5 - 100 ml

## BagMixer® 400

Food, cosmetic, pharmaceutical,  
clinical/medical, agricultural...

Timer : 30 - 210 sec. or ∞, 8 strokes/s.  
40 x 27 x 26 cm, 16.5 kg - 110 V / 220 V  
Capacity : 50 - 400 ml

## BagMixer® 3500 “JumboMix®”

Large samples : textiles, food, filters (residue),  
environment (soil, waste), sterility testing...

Timer : 30 - 210 s. or ∞, 3 - 6 strokes/sec.  
41 x 50 x 38 cm, 33 kg - 110 V / 220 V  
Capacity : 400 - 3500 ml

\*with  
validation  
card

All sterile bags by interscience, BagFilter®, BagPage®, BagLight® fit all models of lab blender on the market.

Please send information to the following address,  
or correct your address, if necessary.

Company \_\_\_\_\_

Name \_\_\_\_\_ First name \_\_\_\_\_


Position \_\_\_\_\_

Address \_\_\_\_\_

Zip code/City \_\_\_\_\_

Country \_\_\_\_\_

Tel. \_\_\_\_\_ Fax \_\_\_\_\_

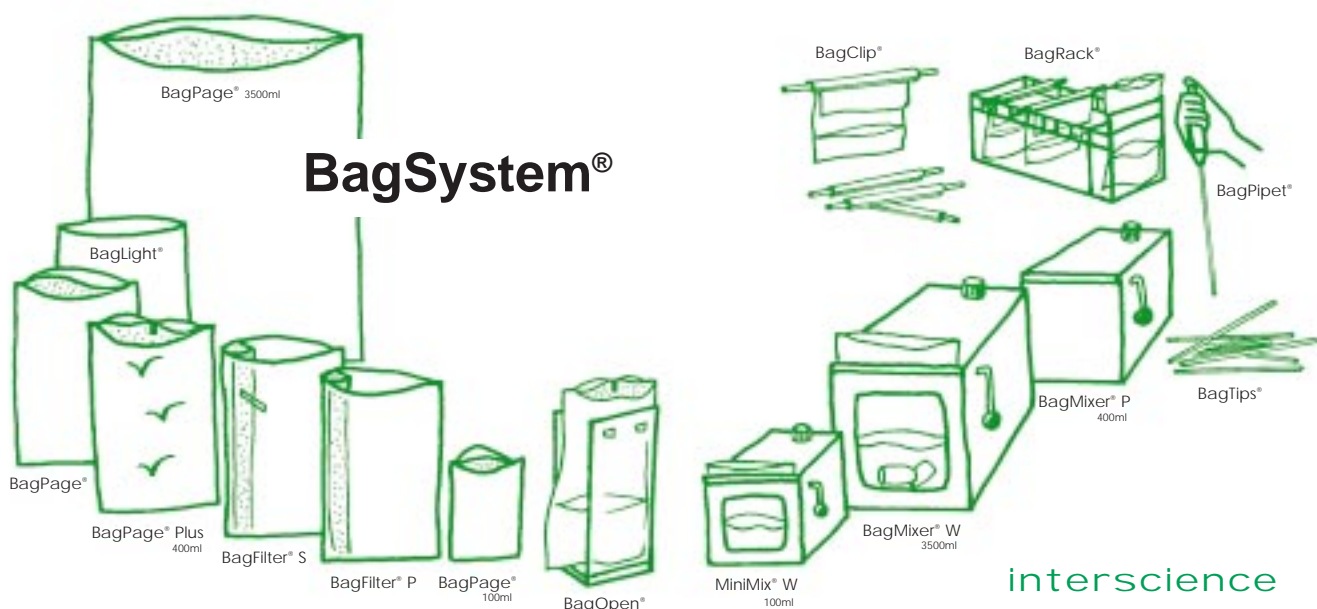

Please circle the products of interest above, or complete this page for more information,  
samples or a demonstration. For a quicker reply, **please fax** this form back to us.



We also propose an attractive range of hygiene products and automats for laboratory. Please consult us.

### Hygiene & Security

	information	samples
<b>Clinet® 101</b> (chlorhexidine disinfecting wipes)	<input type="checkbox"/>	<input type="checkbox"/>
<b>Clinet® 303</b> (disinfecting spray)	<input type="checkbox"/>	<input type="checkbox"/>
<b>Anabac®</b> (autoclave deodorant)	<input type="checkbox"/>	<input type="checkbox"/>

### Automats

	information	demonstration
<b>Spiral® Plater</b> (bacterial enumeration)	<input type="checkbox"/>	<input type="checkbox"/>
<b>Scan® 500</b> (colony counter)	<input type="checkbox"/>	<input type="checkbox"/>
<b>Gravimat®</b> (gravimetric diluter)	<input type="checkbox"/>	<input type="checkbox"/>

your local dealer

KJ 9875 GB

Please return this form to your local dealer or fax to :

**interscience**

100 rte de Ste Gemme  
Tel. +33 (0)1 34 62 62 61  
E-mail : info@interscience.fr

78860 St Nom - France  
Fax +33 (0)1 34 62 43 03  
Web : www.interscience.fr

(Eng.)